

## **510(k) Summary of Safety and Effectiveness**

The following section is included as required by the Safe Medical Device Act (SMDA) of 1990.

Name: pfm Produkte für die Medizin AG  
Address: Wankelstr. 60  
50996 Cologne  
Germany  
CONTACT PERSON: Salvatore F. Palomares, RAC

**JAN 24 2002**

### 510(k) Summary of Safety and Effectiveness

This summary of 510(k) safety and effectiveness information is being submitted in accordance with the requirements of SMDA 1990 and 21 CFR 807.92.

The assigned 510(k) number is: K011783  
Trade Name: Multi-Snare Set  
Common Name: Embolectomy Catheter  
Classification Name: Same

#### Equivalent Devices:

Manufacturer: Microvena Corporation  
Name: Microsnare  
510(k) #: K925439  
Manufacturer: Cook Inc.  
Name: Welter Snare Retriever  
510(k) #: K920823

#### Device Description:

The Multi-Snare Set is designed for the following indications:

- Retrieval and manipulation of foreign objects from the vascular system and hollow viscera
- Assistance in creating loops where cross-over-technique is applied
- Repositioning of indwelling venous catheters

The pfm Multi-Snare Set is contains of the following parts:

- Inner wire which forms the loop;
- Coil mounted on the inner wire in the loop section;
- Shaft body consisting of an outside coil and a core wire; and
- Transition from loop wire to shaft body.

The Multi-Snare Set consists of a highly elastic Nitinol wire with a pre-formed loop on its distal end. Because of the snare's highly elastic design, the snare can be introduced through an insertion sheath. The diameter of the snare can be adjusted depending on the position of the sheath. For better manipulation and torque-ability a torquer is provided.

The introducer sheath is designed for interventional insertion of instruments and systems into the vascular system.

The insertion sheath consists of polyethylene tubing, which is curved at its distal end and has a luer adapter at its proximal end. For better visualization the introducer sheath has an additional distal radiopaque marker.

#### Intended Use:

The Multi-Snare Set (and its components) is intended for:

- Retrieval and manipulation of foreign objects from the vascular system and hollow viscera
- Assistance in creating loops where cross-over technique is applied
- Repositioning of indwelling venous catheters
- Assistance in performing venipuncture to obtain access to central vein

#### Biocompatibility:

The materials used to manufacture the Multi-Snare Set comply with the requirements of ISO 10993-1.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

JAN 24 2002

Mr. Salvatore F. Palomares  
Regulatory Consultant for pfm  
PFM Produkte Fur Die Medizen AG  
C/O Salvatore F. Palomares, RAC  
15 Cherokee Street  
Trabuco Canyon, CA 92679

Re: K011783  
Trade Name: Multi-Snare Ser  
Regulation Number: 21 CFR 870.5150  
Regulation Name: Embolectomy Catheter  
Regulatory Class: II  
Product Code: 74 MMX  
Dated: January 17, 2002  
Received: January 18, 2002

Dear Mr. Palomares:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set

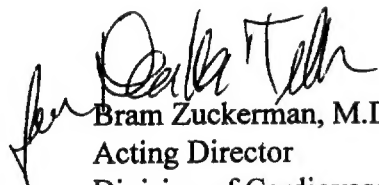
Page 2 - Mr. Salvatore F. Palomares

forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 21 CFR Part 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4586. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read "Bram Zuckerman", is written over the typed name.

Bram Zuckerman, M.D.

Acting Director

Division of Cardiovascular and

Respiratory Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

**510(k):** K011783

**Device Name:** Multi-Snare Set

**Indications for Use:** The Multi-Snare Set (and its components) is intended for:  
Retrieval and manipulation of foreign objects from the vascular system and hollow viscera  
Assistance in creating loops where cross-over technique is applied  
Repositioning of indwelling venous catheters  
Assistance in performing venipuncture to obtain access to central vein

---

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ✓ Or        Over the Counter Use  
(Per 21 CFR 801.109)

  
(Division Sign, On)  
Division of Dental, Infection Control, and  
General Hospital Devices

510(k)  
NUMBER K011783

  
Division of Cardiovascular & Respiratory Devices  
510(k) Number K011783