

JAN 04 2002

## Section J: 510(k) Summary

PORGES™ FOLYSIL™ silicone Foley catheter 510(k) submission

Ref. US1AA61A.DOC

Origin : Regulatory Affairs

K013174  
Page 2 of 2



This 510(k) summary of safety and effectiveness information is being submitted in accordance with the requirements of SDMA 1990.

### J.1. Submitter's information

Submitter's name: PORGES S.A.  
Submitter's address: Centre d'Affaires La Boursidière  
92357 Le Plessis Robinson – France  
Contact person: Mr Bernard ISMAEL  
Regulatory Affairs Manager  
Contact person's phone: + 33 1 46 01 32 06  
Contact person's fax: + 33 1 46 01 32 56  
Contact person's email: [bernard.ismael@porges.com](mailto:bernard.ismael@porges.com)  
Date of preparation: July 2001

### J.2. Device name

Classification name: Urological Catheter, Retention type, Balloon (78 EZL)  
Common / Usual name: Foley catheter  
Proprietary name: PORGES™ FOLYSIL™ silicone Foley catheter

### J.3. Predicate devices

The PORGES™ FOLYSIL™ silicone Foley catheter is substantially equivalent to the modified paediatric silicone Foley catheter from PORGES, silicone Foley catheters from ROCHESTER, Inc., Bardex® All-silicone 3-way Foley catheter from BARD and Oligon® Foley Catheter from IMPEMED.

### J.4. Description of the Device

The PORGES™ FOLYSIL™ silicone Foley catheter is of the retention type, commonly called a Foley catheter. The device is a single use, disposable, sterile with retention balloon, which is attached to the silicone shaft. One lumen is for draining fluids to and from the urinary tract. The second lumen is to inflate and deflate the balloon with sterile water. On models with a third lumen, it is used in conjunction with the first lumen for flushing the urinary tract. Sterile water is used for inflation and deflation of the balloon. The distal end has opposite eye holes, which are used for drainage. On the Tiemann tip there is one hole. Straight cylindrical, "over the guide wire" and Tiemann tips are available versions of this product. On the opposing end of the shaft, are a connecting funnel and a Luer activated valve. This product is available in sizes 12 Fr to 24 Fr, smooth and profile (grooved) shafts. Profile shafts provide drainage channels along the shaft length.

### J.5. Intended use of the Device

The 2-way PORGES™ FOLYSIL™ silicone Foley catheter is used for urethral urinary catheterisation for short term drainage of vesical urines. Only straight 2-way PORGES™ FOLYSIL™ silicone Foley catheter with a maximum 15 ml balloon volume may be used in a supra-pubic way (except the grooved catheters). The 3-way PORGES™ FOLYSIL™ silicone Foley catheter is used for urethral urinary catheterisation to allow short term drainage of the vesical urines, irrigation/injection following surgery.

### J.6. Technological characteristics

The PORGES™ FOLYSIL™ silicone Foley catheter has similar technological and performance characteristics to the predicate devices. The catheter is manufactured entirely from silicone elastomer as for the predicate devices. The catheter is supplied in French sizes ranging from 12 to 24 and balloon capacities 5 cm<sup>3</sup> to 30 cm<sup>3</sup>. The predicate devices are available in French sizes from 6 to 26 balloon capacities from 1.5 cm<sup>3</sup> to 30 cm<sup>3</sup>. The device is supplied in male and female lengths. The predicate devices are supplied in male and female lengths. All of the devices are supplied sterile for single use.

### J.7. Testing and results

The PORGES™ FOLYSIL™ silicone Foley catheter referenced in this submission is held to the same design, manufacture, and performance specifications as those Foley catheters currently manufactured. Performance

K013174  
Page 2 of 2



## Section J: 510(k) Summary

PORGES™ FOLYSIL™ silicone Foley catheter 510(k) submission

Ref. US1AA61A.DOC

Origin : Regulatory Affairs

and functional testing standards are based on the FDA "Guidance for the content of premarket notifications for conventional and antimicrobial Foley catheters" dated September 12, 1994.

The PORGES™ FOLYSIL™ silicone Foley catheter meets the following performance requirements per testing conducted according to ASTM F 623-89, when appropriate, and/or PORGES testing/acceptance criteria:

Note: ASTM F 623-89 excludes from the scope the catheters that have three lumens, balloons equal to or larger than 30 cm<sup>3</sup>, or shaft sizes smaller than 12 Fr or larger. Also excluded are catheters for paediatric and nonurethral catheterizations such as suprapubic cystostomy. However, the test methods described therein will also be utilised to test tri-lumen (3-way) catheters, balloon catheters equal to or larger than 30 cm<sup>3</sup>, and catheters size equal to or smaller than 12 Fr including paediatric catheters.

- Flow rate through the drainage lumen
- Resistance of the balloon to rupture when inflated to the claimed balloon volume and held under conditions approximating the usage environment for a period of seven days;
- Resistance of the inflated balloon to being distorted and pulled through the bladder outlet;
- Maintenance of balloon inflation to fill volume over an extended time;
- Manufacturing tolerances for catheter tip, balloon and shaft diameters;
- Ability of an inflated catheter that has been submerged for seven days to deflate reliably to within 4 Fr. sizes of the labelled shaft size, as applicable, including the time for such deflation;
- Shaft tensile strength and tip adherence;
- Balloon burst.

The PORGES™ FOLYSIL™ silicone Foley catheter passes biocompatibility testing per ISO 10993-1.



Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

JAN 04 2002

Mr. Bernard Ismael  
Regulatory Affairs Manager  
PORGES-C.A. La Boursidière  
92357 Le Plessis Robinson  
CEDEX FRANCE

Re: K013174  
Trade/Device Name: PORGEST™ FOLYSIL™ Silicone  
Foley Catheter  
Regulation Number: 21 CFR 876.5130  
Regulation Name: Urological catheter and accessories  
Regulatory Class: II  
Product Code: 78 EZL  
Dated: December 20, 2001  
Received: December 27, 2001

Dear Mr. Ismael:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (sections 531-542 of the Act); 21 CFR 1000-1050.


This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at one of the following numbers, based on the regulation number at the top of this letter:

8xx.1xxx	(301) 594-4591
876.2xxx, 3xxx, 4xxx, 5xxx	(301) 594-4616
884.2xxx, 3xxx, 4xxx, 5xxx, 6xxx	(301) 594-4616
892.2xxx, 3xxx, 4xxx, 5xxx	(301) 594-4654
Other	(301) 594-4692

Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>.

Sincerely yours,

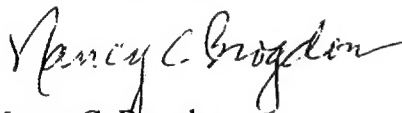

Nancy C. Brogdon  
Director, Division of Reproductive,  
Abdominal, and Radiological Devices  
Office of Device Evaluation  
Center for Devices and Radiological Health

Enclosure

# Section F: Indications for Use Statement

PORGES™ FOLYSIL™ silicone Foley catheter 510(k) submission

Ref. US1AA61A.DOC

Origin : Regulatory Affairs



Page 1 of 1

510(k) Number (if known): K 013174

Device Name: PORGES™ FOLYSIL™ silicone Foley catheter

Indications for use:

The 2-way PORGES™ FOLYSIL™ silicone Foley catheter is used for urethral urinary catheterisation for short term drainage of vesical urines. Only straight 2-way PORGES™ FOLYSIL™ silicone Foley catheter with a maximum 15 ml balloon volume may be used in a supra-pubic way (except the grooved catheters).

The 3-way PORGES™ FOLYSIL™ silicone Foley catheter is used for urethral urinary catheterisation to allow short term drainage of the vesical urines, irrigation/injection following surgery.

(PLEASE DO NOT WRITE BELOW THIS LINE – CONTINUE ON ANOTHER PAGE IF NEEDED)

CONCURRENCE OF CDRH, OFFICE OF DEVICE EVALUATION (ODE)

(Optional Format 3-10-98)

Nancy C. Brydon  
(Division Sign-Off)  
Division of Reproductive, Abdominal,  
and Radiological Devices  
510(k) Number K013174

Prescription Use ✓