

FEB 11 2002

K014301

9.0 SUMMARY OF SAFETY AND EFFECTIVENESS

"510(k) SUMMARY"

- 9.1 Submitter: P/L Biomedical
7690 Cameron Circle
Fort Myers, FL 33912
Tel – 941-768-1118
Fax – 815-550-0162
Contact – Lee Leichter
Prepared – December 12, 2001
- 9.2 Trade/Proprietary Name: Merits Health Products Oxygen Concentrators
- 9.3 Common/Usual Name: Oxygen Concentrator
- 9.4 Classification Name: Portable Oxygen Generator

9.5 Comparison to Currently Marketed Devices

The modified 3-liter Merits Health Products Oxygen Concentrators are substantially equivalent to the currently marketed 5-liter Merits Health Products Oxygen Concentrators (K011884)

9.6 Device Description

The Merits Health Products Oxygen Concentrators are prescription devices designed to provide an inexpensive supply of supplemental oxygen in a home or institution without a continuous source of purified oxygen. They are not life-supporting nor life-sustaining devices. The devices operate through the use of molecular sieve material that binds with the water and nitrogen in filtered room air to leave a gas that is approximately 93% oxygen when delivered to the patient. The compressor creates a vacuum to draw room air into a holding tank. At the same time, downstream of the compressor, the air from the previous cycle is pressurized into one of the two aluminum welded molecular sieve tanks. As the oxygen is forced out of the end of the tank, it enters a "T" fitting that directs most of the gas to flush the nitrogen out of the second molecular sieve tank into the ambient air. The remaining oxygen is delivered to the patient. On the next cycle, the air is directed into the second molecular sieve tank with the oxygen generated flushing the first tank and continuing the supply to the patient. This repetitive cycle generates the oxygen necessary to flush and prepare the saturated sieve tank while supplying the patient with a continuous flow of high concentration oxygen. Options will include an Oxygen alarm and a pediatric flowmeter

9.7 Indications for Use

The oxygen concentrators are intended to provide supplemental oxygen. The device is not intended for life support nor does it provide any patient monitoring capabilities.

9.0 SUMMARY OF SAFETY AND EFFECTIVENESS (Con't)

9.8 Technological Characteristics

The oxygen concentrator operates by using molecular sieve material to absorb water and nitrogen from filtered air. The resulting gas has an increased concentration of oxygen. This technology is well established and has been used in other legally marketed products. There are no major technological differences.

9.9 Performance Data

The results of the testing confirm that the device meets specifications and is substantially equivalent to the predicate device.

9.10 Conclusion

Based on the design, performance specifications and testing and intended use, the Oxygen Concentrators are substantially equivalent to the currently marketed device.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

FEB 11 2002

Mr. Lee Leichter
Merits Health Products Co., Ltd.
c/o P/L Biomedical
7690 Cameron Circle
Fort Myers, FL 33912

Re: K014301
Oxygen Concentrators
Regulation Number: 868.5440
Regulation Name: Generator, Oxygen, Portable
Regulatory Class: II (two)
Product Code: 73 CAW
Dated: January 28, 2002
Received: January 29, 2002

Dear Mr. Leichter:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

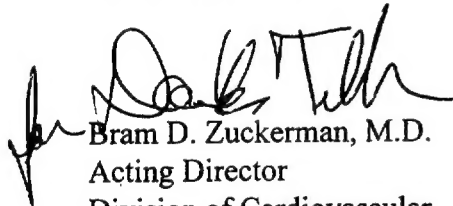
If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 21 CFR Part 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4646. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read "Bram D. Zuckerman".

Bram D. Zuckerman, M.D.
Acting Director
Division of Cardiovascular
and Respiratory Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

INDICATIONS FOR USE STATEMENT

510(k) File Number:

K014301

Device Name:

Merits Health Products Oxygen Concentrator

Indications For Use:

The Oxygen Concentrators are indicated for the delivery of supplemental oxygen in the home or medical institutions. The device is not intended for life support nor does it provide any patient monitoring capabilities.

PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IF NEEDED

Concurrence of CDRH, Office of Device Evaluation (ODE)


Division of Cardiovascular & Respiratory Devices
510(k) Number K014301

Prescription Use ☒
(Per 21 CFR 801.19)

OR

Over-The-Counter Use ☐