



APR 30 2002

1020350

**BCI, Inc.**

N7 W22025 Johnson Road  
Waukesha, WI 53186-1856 USA  
Tel: 262-542-3100  
Fax: 262-542-3325  
www.smiths-bci.com

## **Summary of Safety and Effectiveness**

Submitter:	BCI, Inc.
Address:	N7 W22025 Johnson Road Waukesha, WI 53186
Telephone:	(262) 542-3100
Contact:	VP Regulatory Affairs
Prepared:	January 30, 2002
Proprietary Name:	BCI 3180 Pulse Oximeter
Common/Classification Name:	Pulse Oximeter
Predicate Devices:	BCI 3100 Pulse Oximeter (K893877) BCI 3304 Pulse Oximeter (K962156) BCI 3401 Pulse Oximeter (K980714)

### **New Device Description:**

The BCI 3180 pulse oximeter is a new monitor with the same parameters as existing devices legally marketed by BCI, Inc. This device is designed to provide full featured monitoring capabilities in a tabletop design. The system features an SpO<sub>2</sub> sensor interface, display of patient and waveform data via an LCD display, power status LED, and the function keypad area consisting of six keys (on/off, waveform / trend, alarm audio pause/off, menu / enter, up arrow & down arrow). The monitor has an infrared serial port that is used for data communications to a printer and can download trend data to a printer or personal computer through connection to the SpO<sub>2</sub> sensor interface (sensor must be removed from the device).

**smiths**

Smiths Medical - a part of Smiths Group plc



#### Intended Use:

The BCI 3180 pulse oximeter is indicated for continuous patient monitoring of patient SpO<sub>2</sub>, pulse rate, and pulse strength measurements. It is intended for the hospital or clinical environment and during emergency land transport. It is intended for use in critical environments, including ventilatory applications, patient transport, and anesthesia. It is not intended for use in the home. The pulse oximeter is intended for use with all BCI pulse oximetry sensors, providing SpO<sub>2</sub> and pulse rate measurements on all patients from neonate to adult.

#### Performance Data:

The design of this device utilizes currently available technology found in legally marketed devices. The BCI 3180 Pulse Oximeter is a redesign of the BCI 3100 Pulse Oximeter. The new design uses currently available parts, meet the most recent product standards, and has a more up-to-date case and display. Testing was done to ensure that the BCI 3180 Pulse Oximeter would perform safely and accurately within the environment(s) for which it is to be marketed.

Safety and environmental testing was conducted in accordance with the *Reviewer's Guidance for Respiratory Devices*, 1993, EN 60601-1: 1990, and EN 60601-1-2: 1993. Electromagnetic compatibility (EMC), electrical, mechanical durability, safety (operator and patient), and temperature/humidity testing have been completed. The results demonstrate that the BCI 3180 pulse oximeter is in compliance with the guidelines and standards referenced in the reviewer's guides and that it performs within its specifications and functional requirements.

Testing of device performance included clinical testing of the SpO<sub>2</sub> parameter and overall software validation. The results demonstrated that the BCI 3180 pulse oximeter performed within its specifications.

The testing described above indicates that there is no functional difference between the operation of the BCI 3180 pulse oximeter and the predicate device. Based on these results, it is our determination that the device is safe, effective, and performs as well as the legally marketed predicate device(s).

This summary of 510(k) safety and effectiveness information is being submitted in accordance with the requirements of SMDA 1990 and 21 CFR 807.92.



Respectfully,

A handwritten signature in black ink, reading "Donald Alexander". The signature is fluid and cursive, with a long horizontal line extending from the end of the name.

Donald Alexander  
VP Regulatory Affairs

BCI is a trademark of BCI, Inc. The symbol ® indicates it is registered in the U.S. Patent and Trademark Office and certain other countries.

**smiths**

Smiths Medical - a part of Smiths Group plc



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

APR 30 2002

Mr. Donald Alexander  
BCI, Inc.  
N7 W22025 Johnson Road  
Waukesha, WI 53186

Re: K020350  
BCI 3180 Pulse Oximeter  
Regulation Number: 870.2700  
Regulation Name: Oximeter  
Regulatory Class: II (two)  
Product Code: 74 DQA  
Dated: January 30, 2002  
Received: February 4, 2002

Dear Mr. Alexander:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

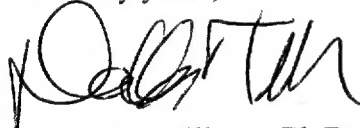
Page 2 – Mr. Donald Alexander

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 21 CFR Part 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4646. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read "Donna-Bea Tillman", written over a horizontal line.

Donna-Bea Tillman, Ph.D.  
Acting Director  
Division of Cardiovascular  
and Respiratory Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

## Indications For Use

510(k) Number (if Known): K020350

Device Name: BCI 3180 Pulse Oximeter

Indications For Use:

Intended Use:

The BCI 3180 pulse oximeter is indicated for continuous patient monitoring of patient SpO<sub>2</sub>, pulse rate, and pulse strength measurements. It is intended for the hospital or clinical environment and during emergency land transport. It is intended for use in critical environments, including ventilatory applications, patient transport, and anesthesia. It is not intended for use in the home. The pulse oximeter is intended for use with all BCI pulse oximetry sensors, providing SpO<sub>2</sub> and pulse rate measurements on all patients from neonate to adult.

( PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IF NEEDED )

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ✓  
(Per 21 CFR 801.109)

OR

Over-The-Counter Use \_\_\_\_\_

  
Division of Cardiovascular & Respiratory Devices  
510(k) Number K020350