



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
10903 New Hampshire Avenue
Document Control Center – WO66-G609
Silver Spring, MD 20993-0002

Ms. Louise M. Focht
Futura Biomedical
990 Park Center Drive, Suite H
Vista, California 92081

December 14, 2016

Re: K033046

Trade/Device Name: Subtalar Peg Implant
Regulation Number: 21 CFR 888.3040
Regulation Name: Smooth or Threaded Metallic Bone Fixation Fastener
Regulatory Class: Class II
Product Code: MJW
Dated: September 29, 2003
Received: September 29, 2003

Dear Ms. Focht:

This letter corrects our substantially equivalent letter of December 23, 2003.

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration. Please note: CDRH does not evaluate information related to contract liability warranties. We remind you, however, that device labeling must be truthful and not misleading.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Parts 801 and 809); medical device reporting (reporting of medical device-related adverse events) (21 CFR 803); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the

electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

If you desire specific advice for your device on our labeling regulation (21 CFR Parts 801 and 809]), please contact the Division of Industry and Consumer Education at its toll-free number (800) 638-2041 or (301) 796-7100 or at its Internet address <http://www.fda.gov/MedicalDevices/ResourcesforYou/Industry/default.htm>. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). For questions regarding the reporting of adverse events under the MDR regulation (21 CFR Part 803), please go to <http://www.fda.gov/MedicalDevices/Safety/ReportaProblem/default.htm> for the CDRH's Office of Surveillance and Biometrics/Division of Postmarket Surveillance.

You may obtain other general information on your responsibilities under the Act from the Division of Industry and Consumer Education at its toll-free number (800) 638-2041 or (301) 796-7100 or at its Internet address <http://www.fda.gov/MedicalDevices/ResourcesforYou/Industry/default.htm>.

Sincerely yours,

Mark N. Melkerson -S

Mark N. Melkerson
Director
Division of Orthopedic Devices
Office of Device Evaluation
Center for Devices and Radiological Health

Enclosure

510 (k) Number (If Known): K033046
Device Name: Subtalar Peg Implant

Page 1 of 1

Indications for Use:

The Futura Biomedical Subtalar Peg Implant is intended to treat hyperpronated foot and stabilize the subtalar joint. It is intended to block forward, downward and medial displacement of the talus, thus allowing normal subtalar joint motion but blocking excessive pronation and the resulting sequela. Examples include:

- Flat foot treatment in children and adolescents
- Congenital flat foot
- Unsuccessful long term orthopaedic treatment (shoes, insoles...)
- Tarsal coalitions
- Painful flat foot
- Supple deformity in posterior tibial tendon dysfunction
- Paralytic flat foot
- Subtalar instability

This device is intended to be fixed with bone cement.


Prescription Use X
(Part 21 CFR 801 Subpart D)

AND/OR

Over-The-Counter Use _____
(21 CFR 807 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

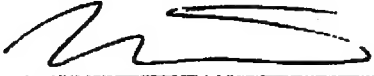

Deputy Associate Director
Office of Device Evaluation

K033046

DEC 23 2003

K033046

PAGE 1 of 3

510 (k) Summary

Summary of 510 (k) safety and effectiveness information upon which the substantial equivalence determination is based:

Prepared: December 22, 2003
Applicant: Futura Biomedical
990 Park Center Drive
Vista, CA 92081
Telephone: 760-599-1670
Fax: 760-599-1675
Contact: Louise M. Focht

| | |
|------------------------|-----------------------------|
| Device Name: | Subtalar Peg Implant |
| Device Trade Name: | Subtalar Peg Implant |
| Device Classification: | Class II |
| Reviewing Panel: | Orthopedic |
| Regulation Number | |
| Product Code: | 87 MJW |
| Predicate Device: | K792670- Smith Subtalar Peg |
| Registration Number: | 2030833 |
| Owner Operator Number: | 9028319 |

Device Description:

The Futura Biomedical Implant is a one-piece device made of UHMWPE or Cobalt Chromium with titanium plasma spray on the stem, intended to be implanted into the calcaneus of the foot. The implant is designed in 13 sizes and two configurations. The implant which is used in the treatment of excessive motion of the talus relative to the calcaneus acts as a spacer for the joint, maintaining the joint space, allowing for range of motion, but limiting excessive pronation.

Indications for Use:

The Futura Biomedical Subtalar Peg Implant is intended to treat hyperpronated foot and stabilize the subtalar joint. It is intended to block forward, downward and medial displacement of the talus, thus allowing normal subtalar joint motion but blocking excessive pronation and the resulting sequela. Examples include:

- Flat foot treatment in children and adolescents
- Congenital flat foot
- Unsuccessful long term orthopaedic treatment (shoes, insoles...)
- Tarsal coalitions

Painful flat foot
Supple deformity in posterior tibial tendon dysfunction
Paralytic flat foot
Subtalar instability

This device is intended to be fixed with bone cement.

Comparison to Predicate Device:

The legally marketed predicate device to which this device is substantially equivalent is the Wright Medical STA-Peg implant.

Regulatory Class: II
Product Code: 87 MJW

Table 1. Comparison of Futura Biomedical and Wright Medical STA-Peg Implant

| Item | Futura Product | Wright Medical Product |
|---------------------|---|---|
| Product Name | Subtalar Peg Implant | STA-Peg Implant |
| Use | Single use | Single use |
| Fixation | Stem in bone | Stem in bone |
| Constraint | Non constrained | non constrained |
| Material | UHMWPE or CoCr/CpTi | UHMWPE |
| Sizes | 5 sizes, two configurations (13 sizes) | 5 sizes |
| Indications for use | <p>The Futura Biomedical Subtalar Peg Implant is intended to treat hyperpronated foot and stabilize the subtalar joint. It is intended to block forward, downward and medial displacement of the talus, thus allowing normal subtalar joint motion but blocking excessive pronation and the resulting sequela. Examples include:</p> <p>Flat foot treatment in children and adolescents Congenital flat foot Unsuccessful long term orthopaedic treatment (shoes, insoles...) Tarsal coalitions Painful flat foot Supple deformity in posterior tibial tendon dysfunction Paralytic flat foot Subtalar instability</p> <p>This device is intended to be fixed with bone cement.</p> | <p>The STA-Peg implant is intended to treat hyperpronated foot and stabilize the subtalar joint. It is intended to block forward, downward and medial displacement of the talus, thus allowing normal subtalar joint motion but blocking excessive pronation and the resulting sequela. Examples include:</p> <p>Flat foot treatment in children and adolescents Congenital flat foot Unsuccessful long term orthopaedic treatment (shoes, insoles...) Tarsal coalitions Painful flat foot Supple deformity in posterior tibial tendon dysfunction Paralytic flat foot Subtalar instability</p> |

Similarities of the Futura Biomedical Subtalar Arthrorisis Implant and the Wright Medical STA-Peg Implant include:

Both devices are: intended for single use only; intended for surgical implantation longer than 30 days; both devices are placed into the calcaneus of the foot, allowing normal subtalar joint motion while blocking excessive pronation and resulting sequela; both devices are made of industry standard materials, no new materials are introduced in either product; Both devices are comparably sized; both devices have the same indications for use.

Summary:

The device and the predicate device have similar design characteristics and intended use. The new device is substantially equivalent to the predicate device.